


# LET’S MOVE AGAIN.
# OUTDOOR ROLLATOR.

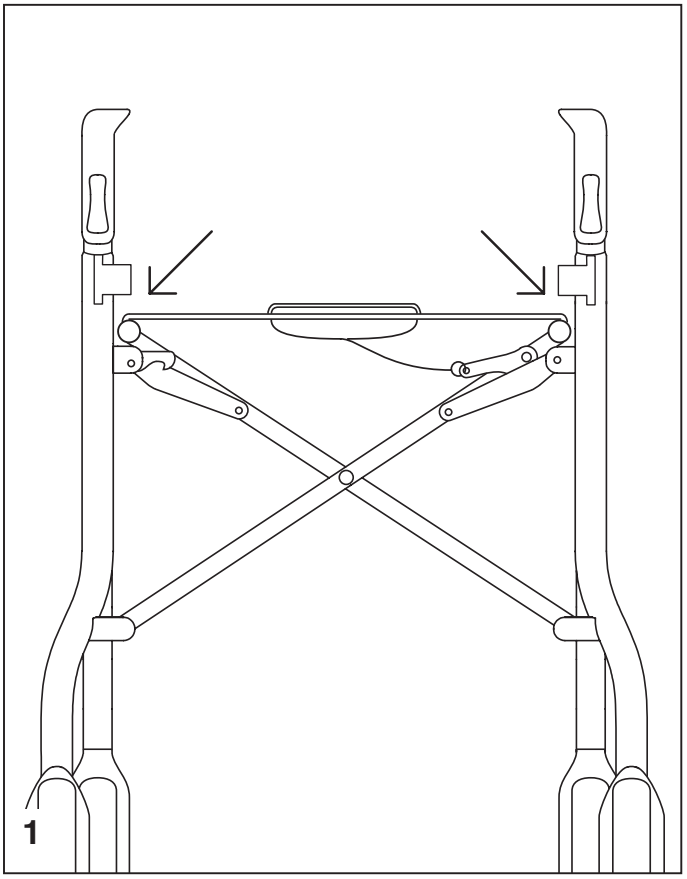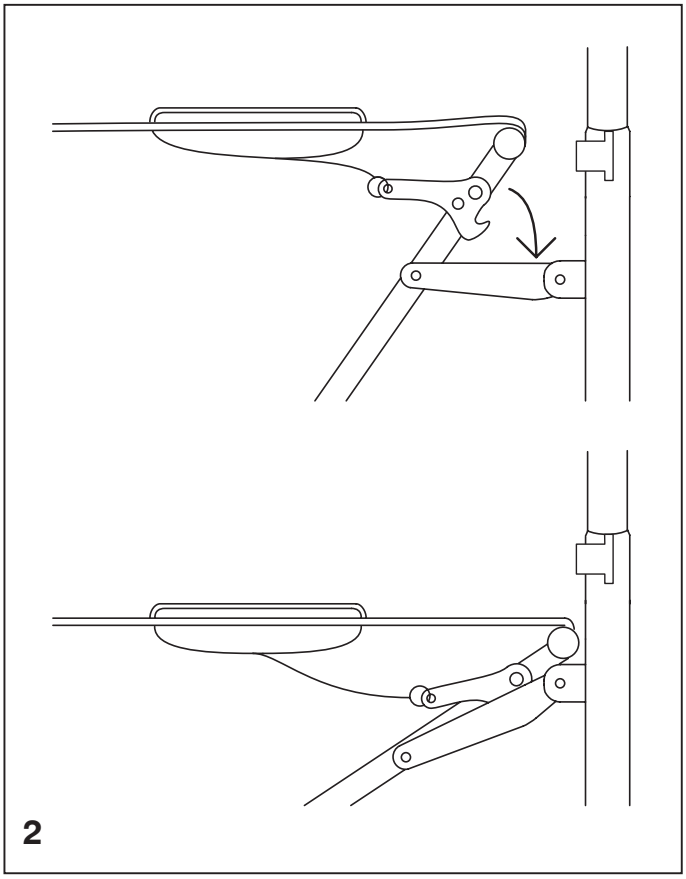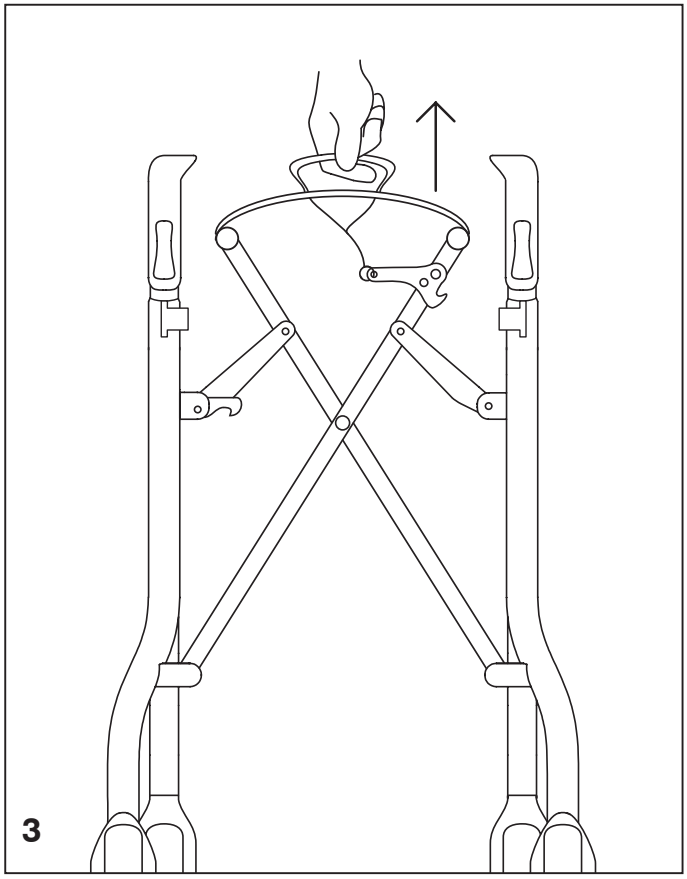
User Manual
June 2026

1


2

3

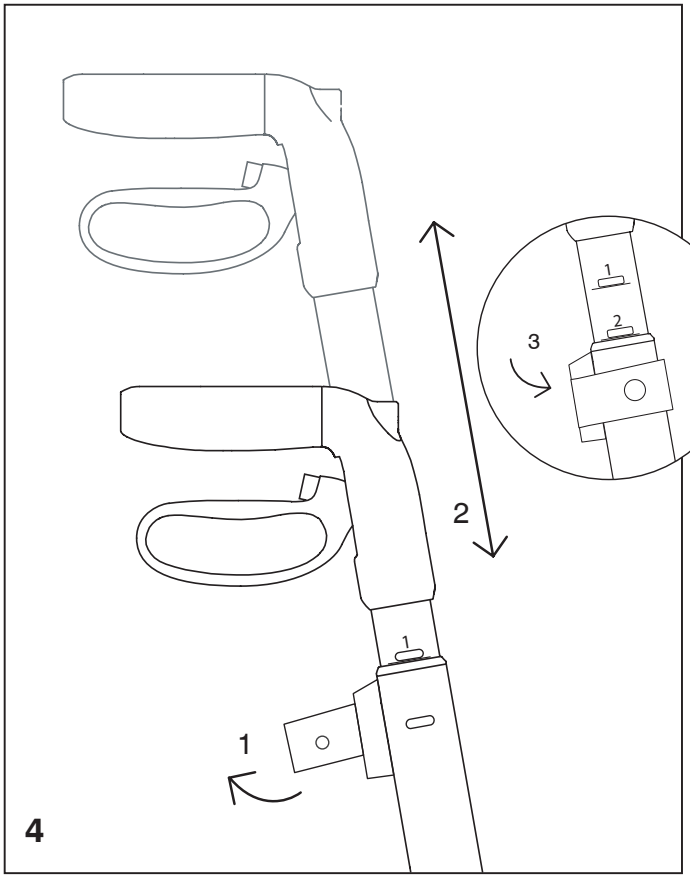1
2
3
1
2
1
4

## MATERIAL

01. Handgrip - TPR
02. Handle Tube - Aluminum
03. Seat Holders - Steel
04. Frame - Aluminum
05. Front Fork - PA6 + 30% GF
06. Rear Fork - PA6 + 30% GF
07. Tires - TPR
08. Wheel rim - PP
09. Seat - PVC + Nylon
10. Bushing for handle - PA
11. Folding cross - Steel

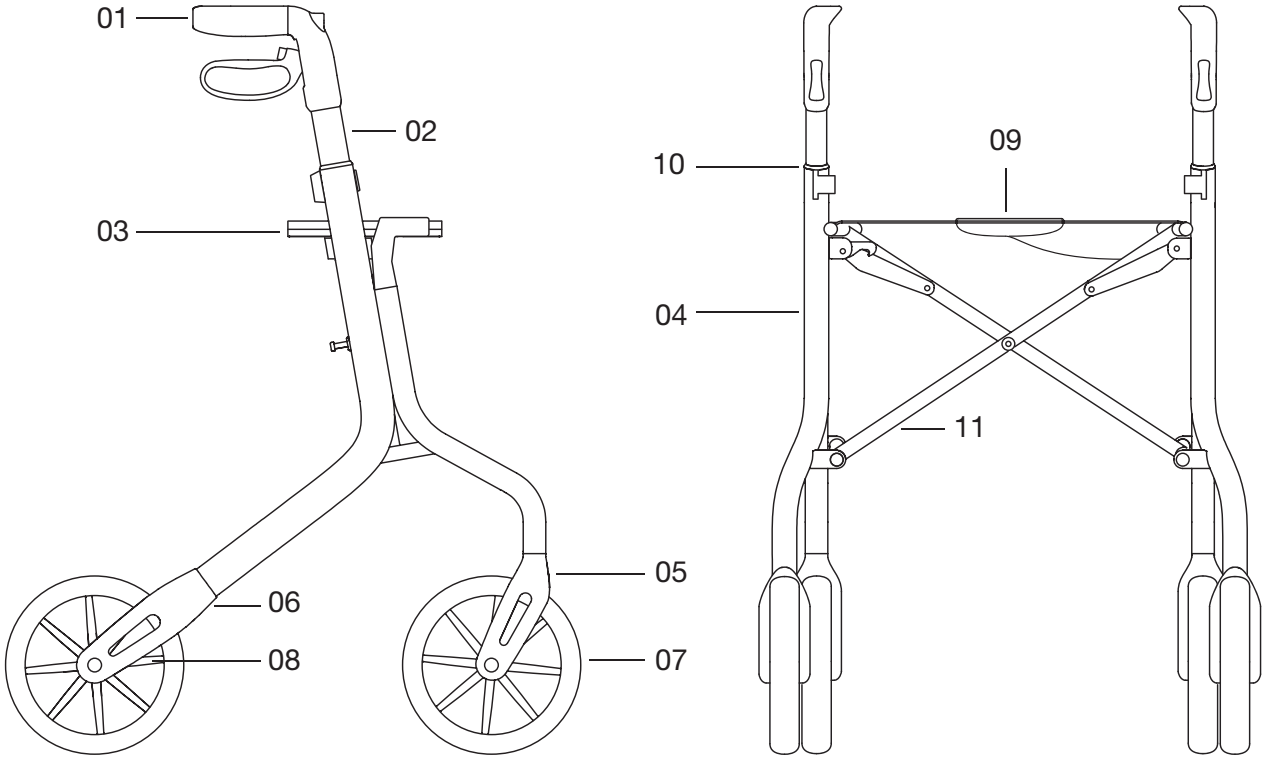**Återvinning / Genbrug / Recycling / Recycling / Recycling / Riciclaggio / Recyclage**

Denna produkt skall återvinnas i enlighet med nationella riktlinjer.
Dette produkt skal genbruges i overensstemmelse med nationale retningslinjer.
This product should be recycled in accordance with national guidelines.
Dit product moet worden gerecycled in overeenstemming met de nationale richtlijnen. Dieses Produkt sollte gemäß den nationalen Richtlinien recycelt werden.
Questo prodotto deve essere riciclato in conformità con le linee guida nazionali.
Ce produit doit être recyclé conformément aux directives nationales.

# ENGLISH

**Let's Move Again Outdoor Walker.**

We congratulate you on your new rollator, Let's Move Again, that will make your life a little easier, more pleasant and more safe. To make it easier for you to get familiar with your rollator it is important that you study this manual. The rollator helps you to walk and move outdoors. You can also use it to sit on when you want to rest during a walk. Then do not forget to activate the parking brake – position 3 on the picture to the right. To cross curbs and other obstacles, walk in a straight angle to the curb for the rollator not to become unsteady as it crosses the edge. The seat is designed only to sit on, not for transport of for instance children or objects. Do not put heavy objects like large boxes, pieces of furniture etc. on the seat.

Fold the rollator for transport. Ensure that the rollator is properly secured during transport. Take extra precautions to secure the equipment during transport, if necessary.


**Maintenance**
Clean the rollator with lukewarm water and a little detergent. Don't put the rollator in hot places or where it is very damp.

**Warranty**
We provide a two-year warranty. The warranty covers manufacturing defects and does not apply in cases of abnormal wear or misuse. Wear parts such as cables and wheels are not covered by the warranty.


**How to unfold**
(See pictures on page 2!)
1. Rollator standing. Unfold the rollator by pushing the two sides apart from each other.
2. Make sure that the rollator is fully unfolded and the locking hook under the seat snaps in place.

**How to fold**
3. Lift the seat in the middle and the locking hook under the seat is released and the rollator is folded. Press the sides together if needed.

**Handle height adjustment**
4. Release the adjustment lock, align the height with the tube's indicators, and fasten the lock back. (Picture 4 page 2!)

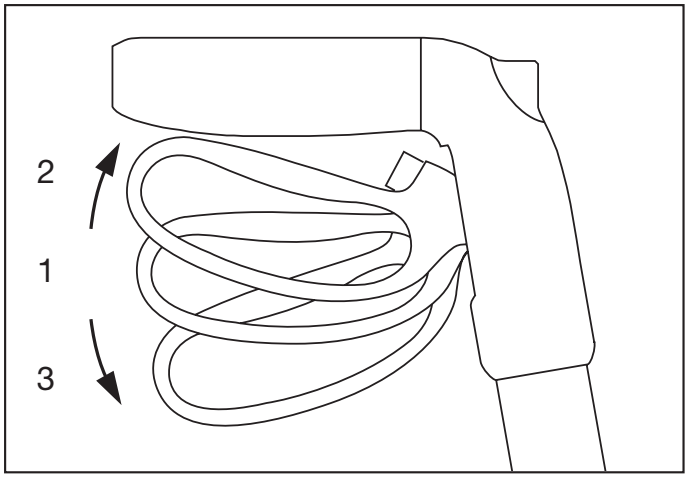**Brake functions**
1. In position 1 the rollator is rolling freely.
2. To brake the rollator pull the handle up to position 2.
3. For parking brake push handle down to position 3, where the brake locks.

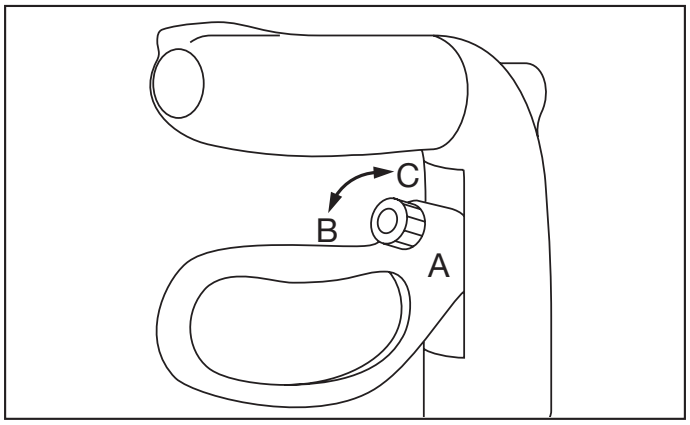**Brake adjustment**
1. The brakes are adjusted with the adjustment screw at the brake handle (A).
2. Adjustment of the wire tension is done by turning the adjustment screw (B-C).
If the brake is too loose, turn in direction (C). If the brake is to hard turn in direction (B).

Tested according to ISO 11199-2:2005
for a maximum user weight of 130 kg.

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